

EP14

Evaluation of Commutability of Processed Samples

Sample

This document provides guidance for evaluating the commutability of processed samples by determining if they behave differently than unprocessed patient samples when two quantitative measurement procedures are compared.

A guideline for global application developed through the Clinical and Laboratory Standards Institute consensus process.

Evaluation of Commutability of Processed Samples

Karl De Vore, BA, SSBB
Yashpal Agrawal, MD, PhD
Todd D. Alspach, MT(ASCP), MOS
Jeffrey R. Budd, PhD
Ramon A. Durazo-Arvizu, PhD
John H. Eckfeldt, MD, PhD
Kathie Goodwin, RAC, MBA, MT(ASCP)BB

Abdel-Baset Halim, PharmD, PhD, DABCC
Thomas A. Long
W. Gregory Miller, PhD
Nisar Pampori, PhD
Justin Thomas
Jeffrey E. Vaks, PhD
Hubert W. Vesper, PhD

Abstract

Clinical and Laboratory Standards Institute document EP14—*Evaluation of Commutability of Processed Samples* was developed for manufacturers and providers of proficiency testing or external quality assessment programs, although it is useful to clinical laboratories as well. The document helps users 1) determine whether noncommutability is the source of unexpected results that are sometimes observed with processed samples when two quantitative measurement procedures are compared, 2) display the magnitude of the effects, and 3) ensure that laboratory performance is evaluated fairly if noncommutability is present. The suggested protocol was developed using unprocessed patient samples as the standard of comparison.

Clinical and Laboratory Standards Institute (CLSI). *Evaluation of Commutability of Processed Samples*. 4th ed. CLSI guideline EP14 (ISBN 978-1-68440-156-7 [Print]; ISBN 978-1-68440-157-4 [Electronic]). Clinical and Laboratory Standards Institute, USA, 2022.

The Clinical and Laboratory Standards Institute consensus process, which is the mechanism for moving a document through two or more levels of review by the health care community, is an ongoing process. Users should expect revised editions of any given document. Because rapid changes in technology may affect the procedures, methods, and protocols in a standard or guideline, users should replace outdated editions with the current editions of CLSI documents. Current editions are listed in the CLSI catalog and posted on our website at www.clsi.org.

If you or your organization is not a member and would like to become one, or to request a copy of the catalog, contact us at:

P: +1.610.688.0100 **F:** +1.610.688.0700 **E:** customerservice@clsi.org **W:** www.clsi.org

Copyright ©2022 Clinical and Laboratory Standards Institute. Except as stated below, any reproduction of content from a CLSI copyrighted standard, guideline, companion product, or other material requires express written consent from CLSI. All rights reserved. Interested parties may send permission requests to permissions@clsi.org.

CLSI hereby grants permission to each individual member or purchaser to make a single reproduction of this publication for use in its laboratory procedure manual at a single site. To request permission to use this publication in any other manner, e-mail permissions@clsi.org.

Suggested Citation

CLSI. *Evaluation of Commutability of Processed Samples*. 4th ed. CLSI guideline EP14. Clinical and Laboratory Standards Institute; 2022.

Previous Editions:

Proposed Guideline

April 1999

Approved Guideline

March 2001

Approved Guideline—Second Edition

January 2005

Approved Guideline—Third Edition

August 2014

EP14-Ed4

ISBN 978-1-68440-156-7 (Print)

ISBN 978-1-68440-157-4 (Electronic)

ISSN 1558-6502 (Print)

ISSN 2162-2914 (Electronic)

Volume 42, Number 17

Contents

Abstract	i
Committee Membership	iii
Foreword	vii
Chapter 1: Introduction	1
1.1 Scope	2
1.2 Background	3
1.3 Standard Precautions	5
1.4 Terminology	5
Chapter 2: Commutability Assessment	9
2.1 Assessment Methods	10
2.2 Process Flow Chart	12
2.3 Materials and Samples Assembly	13
2.4 Measurement Procedures	14
2.5 Data Analysis	14
2.6 Documenting Results of the Commutability Evaluation	20
Chapter 3: Conclusion	21
Chapter 4: Supplemental Information	23
References	24
Additional Resources	26
Appendix A. Description of Mathematical Model Used for Evaluating Commutability of Processed Samples Using Deming Regression	27
Appendix B. Outlier Evaluation for a Measurement Procedure Comparison Using Deming Regression	33
Appendix C. Examples of Completed Analyses	35
The Quality Management System Approach	51

Foreword

The primary goal of manufacturers of *in vitro* diagnostic (IVD) measurement procedures is to accurately report measurand values for patient samples. These measurement procedures may not produce accurate results when processed samples such as external quality assessment (EQA) samples, proficiency testing (PT) samples, or QC samples are measured. Because such processed sample matrixes typically undergo supplementation with additional components, and therefore are altered in some manner, measured results may not reflect the accuracy that would be observed for patient samples. Processed samples that have the same measurement response as that of patient samples that contain the same amount of an analyte are called commutable, while those that do not are called noncommutable. In this document, a matrix effect is defined broadly as differing test result biases in processed samples vs unprocessed patient samples due to unknown causes. The matrix effects that cause biases compared with unprocessed patient samples could be correlated to differences in conditions as encompassing as the entire measurement system or as specific as a reagent lot within a single IVD device.

Biases due to matrix effects in processed samples have the potential to affect the quality of patient care because of an inaccurate estimation of the accuracy of a measurement procedure. Depending on the intended use of the processed sample, the impact can range from negligible to serious. For example, a specific bias in a measuring interval verification sample set may have a different impact on the quality of patient care than the same bias in a QC sample. A measuring interval sample set matrix-related bias can directly affect the measuring interval allowed in patient sample results, whereas a QC matrix-related bias may affect the interpretation of QC results following a reagent lot change.

The objective of EP14 is to provide methods for identifying noncommutability so that improvements in measurement procedure analytical specificity and matrix compatibility may be considered. For example, a beneficial outcome of an evaluation may be a change in the processed sample's matrix or its additives, with an improvement in sample commutability. The techniques described in this guideline are also valid for testing the commutability of other samples, such as patient samples that have been processed (eg, added preservatives or spiking material, diluted, depleted, or frozen). Such samples, designed to substitute for unprocessed patient samples, are referred to as surrogate samples in CLSI document EP39.¹ EP14 will be helpful in exploring differences in test material results between measurement procedures, especially when such material serves as a basis for determining measurement procedure performance.

Overview of Changes

This guideline was revised in 2022 under the Limited Revision Process and replaces the third edition of the guideline, which was published in 2014. Several changes were made in this edition, including:

- Adding content focused on processed samples such as EQA samples, PT samples, QC samples, and altered patient samples
- Adding an option for using preset commutability assessment criteria instead of prediction intervals
- Clarifying that EP14 is not designed for use with IVD manufacturer-specific calibrators and that the assessment methods described in this guideline should not be used in product regulatory submissions from such manufacturers

KEY WORDS

bias

commutability

Deming regression

interference

matrix

matrix effect

Evaluation of Commutability of Processed Samples

1 Introduction

This section includes:

- Document scope and applicable exclusions
- Background information pertinent to the document content
- Standard Precautions information, as applicable
- Terms and definitions used in the document
- “Note on Terminology” that highlights particular use and/or variation in use of terms and/or definitions, where applicable
- Abbreviations and acronyms used in the document


1.1 Scope

This guideline provides protocols for evaluating commutability of processed samples when tested with quantitative measurement procedures. Such processed samples may be those created for proficiency testing/external quality assessment (PT/EQA), measuring interval verification sample sets, or QC samples. Processed samples can also be human specimens that are modified in a way that may change their measurement characteristics. In such cases, only a few processed samples (eg, three to five) should be evaluated to represent the behavior of the modification process being assessed. This restriction is suggested because the underlying analyses in this guideline do not account for the simultaneous assessment of many processed samples. CLSI document EP39¹ mentions the use of surrogate samples for measurement procedure performance testing, sometimes using many such samples, eg, for measurement procedure comparison studies. EP14 does not provide methods to screen a large number of samples in such cases.

This guideline is intended for developers of commercial diagnostic tests as well as laboratory-developed tests. This guideline is also useful for manufacturers of measuring interval sample sets and QC samples, and for PT or EQA providers to determine whether their processed samples are commutable with patient samples when examined with specified measurement procedures. Other options for assessing the commutability of processed samples across multiple measurement procedures are covered in CLSI document EP30² and in the literature.³

This guideline may also be useful to all clinical laboratory professionals wishing to investigate the commutability of a processed sample analyzed with a specific *in vitro* diagnostic (IVD) device in their laboratory.

EP14 is intended to assist in the education of clinical laboratorians and diagnostic manufacturers about the commutability of processed materials and how a sample’s matrix can affect some measurand values and their interpretation (referred to as matrix effects). For example, professionals may not be warned of a matrix effect caused by the interaction of processed PT/EQA material and the measurement procedure, and therefore the data may suggest to them that erroneous patient results are being generated, when in fact the results may be acceptable. Examples of a matrix effect due to the interaction of a processed QC and certain reagent lot(s) exist in the literature.⁴ Therefore, these types of effects should not be a surprise to experienced laboratory staff and should not lead to erroneous conclusions about the suitability of results for patient samples. This guideline should assist all interested parties in not only evaluating the presence or absence of a matrix effect, but also by increasing awareness that the intended use of a processed matrix potentially affects patient care.

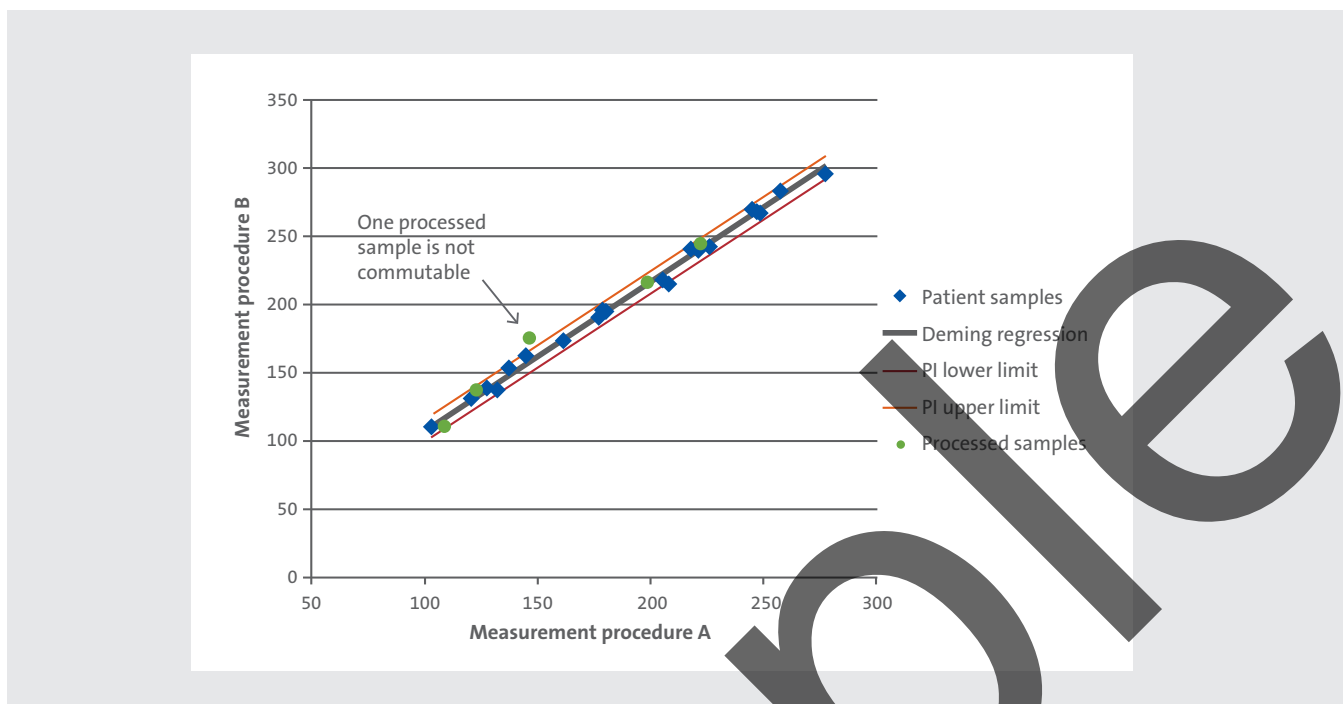

Abbreviation: PI, prediction interval.

Figure 6. Deming Regression of Measurement Procedure B vs Measurement Procedure A

2.5.4 Comparison of Processed Samples With Patient Samples Using Preset Criteria

If preset criteria are used to compare processed samples with patient samples, the comparison methods used should be the same as those described above in Section 2.5.3. For relationships in which the scatter is consistent over the measured interval, a constant offset (in units) above and below the Deming regression line should be used as the preset criterion. That preset criteria may be based on a clinically significant difference. The resultant plot will look similar to Figure 6, with the preset criterion replacing the upper and lower limits of the PI. For relationships in which the scatter increased with increased concentration or activity, a proportional offset (ie, percent difference) above and below the Deming regression line should be used. In this case, the data, including the bias criterion, should be plotted using a log transformation. The resultant plot will look similar to Figure 5, with the preset criterion replacing the upper and lower limits of the PI. Examples of datasets for which both PI and preset limits are applied are provided in Appendix C.

In the ideal situation, the two measurement procedures used in the commutability assessment will have low imprecision and high selectivity for the measurand resulting in a narrow distribution of patient sample points around the regression. Also ideally, there will be an established definition of a difference that is clinically significant. In this case, the preset criteria based on this difference can be used. However, if the measurement procedures used for the commutability assessment have high imprecision and low selectivity, or if a clinically significant difference has not been established, the best option may be to use the PI for the commutability assessment.

Appendix A. Description of Mathematical Model Used for Evaluating Commutability of Processed Samples Using Deming Regression

When a set of patient samples is measured by two different measurement procedures, it is expected that the relationship between the measured values will follow a linear relationship (see CLSI document EP06¹). Under the preceding conditions, the relationships between the results for the two methods can be presented for the patient samples and the processed sample as:

$$Y = \alpha_H + \beta_H X \quad (1a)$$

$$Y = \alpha_{Pc} + \beta_{Pc} X \quad (1b)$$

where:

α_{Pc}, α_H = intercepts,

β_{Pc}, β_H = slopes,

and the subscripts H and Pc indicate patient samples (human) and processed samples, respectively.

Note that the number of replicates in the equations that follow use the single annotation N_H . This implies that the same number of replicates is used for both human (H) and processed samples (Pc). This also implies that the same number of samples is used for both measurement procedures. As described in the document text, if the repeatabilities of the measurement procedures or the sample types are different and this difference is known, then the user may choose to have different numbers of replicates for each situation. In such a case, the averages computed in Equations 4, 6 and 12 should use the number of replicates for each sample determination, rather than a common annotation N_H .

The equivalence of the mathematical relationships in Equations (1a) and (1b) would be established by showing that the respective model parameters are equal pairwise:

$$\alpha_{Pc} = \alpha_H \quad (2a)$$

$$\beta_{Pc} = \beta_H \quad (2b)$$

Equations (2a) and (2b) assume no measurement error. With measurement errors in both measurement procedures, Equation (1a) for the patient samples can be expressed as:

$$Y = \alpha_H + \beta_H (X + \varepsilon_X) + \varepsilon_Y \quad (3)$$

where $\varepsilon_X, \varepsilon_Y$ are random errors in the X and Y measurement procedures. Equation (3) parameters can be estimated with data using regular Deming regression under the following assumptions: the random errors $\varepsilon_X, \varepsilon_Y$ are independent (across the measurement procedures, samples, and replicates) and normally distributed with 0 means and have constant, measurand level-independent SDs, $\sigma(\varepsilon_X), \sigma(\varepsilon_Y)$. If these assumptions do not apply to the obtained data, transformations must be performed as described in Section 2.5 of this document.

Discover How CLSI Can Improve Your Organization



The leading source for the latest medical laboratory standards.



CLSI membership lets you directly impact best practice standards used to improve patient care worldwide—standards you use every day. Membership provides you with standards access, volunteering opportunities, influence in the standards development process, networking opportunities, discounts, and more.

Discover the membership option for you at clsi.org/join.



Our educational and training programs provide convenient, cost-effective continuing education and training resources to help you advance your professional development. We have a variety of easy-to-use, online educational resources and in-person trainings that make learning stress-free and convenient for you and your staff.

See our current offerings at clsi.org/global-training.



Ensure high-quality laboratory testing with CLSI standards. eCLIPSE Ultimate Access™, our complete online library of standards, makes it easy for you and your staff to quickly find the CLSI resources you need. Read, search, link, annotate, bookmark, and share notes with your staff, all within one easy-to-use platform.

Learn more at clsi.org/eCLIPSE.

Sample



CLINICAL AND
LABORATORY
STANDARDS
INSTITUTE®

PRINT ISBN 978-1-68440-156-7

ELECTRONIC ISBN 978-1-68440-157-4

EP14-Ed4